




TAURUS WELLNESS

Operating Instructions

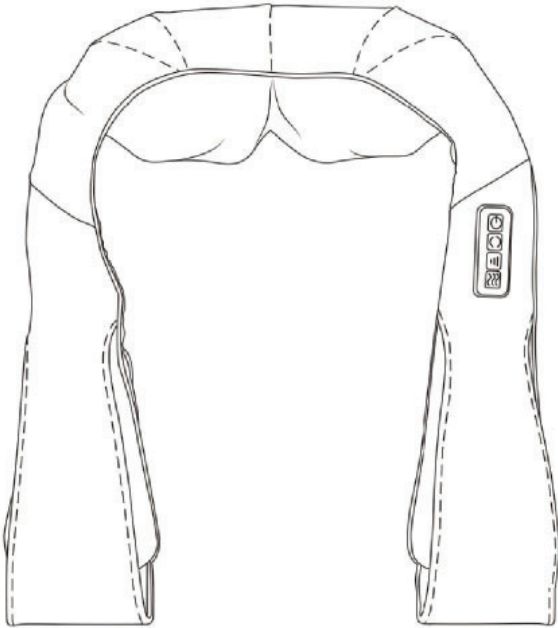

Art. No. TW-NECK

FSTWNECK.01.01

GENERAL SAFETY INSTRUCTIONS

Personal Safety

⚠ DANGER

- + Before you start using the equipment, you should consult your physician that this type of exercise is suitable for you from a health perspective.
- + Note that excessive training can seriously endanger your health. If you notice any signs of weakness, nausea, dizziness, pain, shortness of breath, or other abnormal symptoms, stop exercising immediately and seek advice from your doctor if necessary.

⚠ WARNING

- + This equipment may not be used by children under 14 years.
- + Children should not be allowed unsupervised access to the equipment
- + Persons with disabilities must have a medical license and must be under strict observation when using the equipment.
- + The equipment is strictly for use by one person at a time.

⚠ CAUTION

- + Make sure that nobody is within the range of motion of the equipment during training so as not to endanger you or other persons.

▶ ATTENTION

- + Do not insert any objects of any kind into the openings of the device.

Electrical Safety

⚠ DANGER

- + To reduce the risk of an electric shock, always unplug the equipment from the mains socket before use, maintenance or cleaning. Do not pull on the cable.
- + To reduce the risk of an electric shock, do not use the device if you have previously been in contact with water and/or are in a humid environment.

⚠ WARNING

- + Do not leave the equipment unattended while the mains cable is plugged into the mains socket. During your absence, the mains cable must be removed from the mains socket to prevent improper use by third parties or children.
- + If the mains cable or plug is damaged or defective, contact your contract partner. Until repair, the equipment must not be used.

▶ ATTENTION

- + The equipment requires a mains connection of 220-230 V with 50 Hz mains voltage.
- + The equipment may only be connected directly to an earthed socket using the supplied mains cable. Extension cables must conform to VDE guidelines. Always completely unwind the mains cable.
- + The socket must be protected by a fuse with a minimum fuse rating of „16 A, slow blow“.
- + Do not make any changes to the mains cable or the mains plug.
- + Keep the mains cable away from water, heat, oil and sharp edges. Do not route the mains cable underneath the equipment or under a carpet or rug, and do not place any objects on top of it.
- + Do not unplug the mains cable until the product is turned off.

Set-Up Place

⚠ WARNING

- + Do not place the equipment in main corridors or escape routes.
- + To prevent fire, do not place the device near or in areas with high heat (e.g. near a heater).

⚠ CAUTION

- + Choose the set-up place so that there is sufficient space around the product.

▶ ATTENTION

- + The device may only be used in one building, in sufficiently tempered and dry rooms (ambient temperatures between 10°C and 35°C). The equipment should not be used outdoors or in rooms with high humidity (over 70%) like swimming pools.

Intended Use

The equipment may only be used for its intended purpose.

The equipment is only suitable for home use. The equipment is not suitable for semi-professional (e.g. hospitals, associations, hotels, schools, etc.) and commercial or professional use (e.g. fitness studios).

Technical data

Rated voltage:	12.6V 2A
Rated power:	25.2 W
Li-on battery:	11.1V 2000mAh
Massage duration:	15 minutes
Article weight (gross, incl. packaging):	2.1 kg
Article weight (net, without packaging):	2 kg
Packaging dimensions (L x W x H):	41.5 cm x 19.5 cm x 23.5 cm
Set-up dimensions (L x W x H):	40.5 cm x 18.5 cm x 22.5 cm
Use class:	H



Power button: On/Off



Rotate direction: Clockwise / anticlockwise



Massage speed: High / medium / low



Heating button: On/Off

Neck Massager

Operating Instructions

Charging the Battery

NOTICE

- + The device must be fully charged before first use.
- + Fully charged, the device can be used for 6 massages in a row (15 minutes per massage).
- + The massager cannot be used while charging.

1. Connect the power adapter to the device and then plug it into a power outlet.
2. Disconnect the device from the mains adapter and from the socket when the charging process is finished.

Start Massage

1. To turn on the device and start the massage, press and hold the power button for 3 seconds.

The massage begins at high speed and with the heating on.


The heating function, massage speed and massage direction can be adjusted by pressing the respective buttons.

2. To end the massage and turn off the device, press and hold the power button for 3 seconds.


Massage Areas



Neck & shoulder

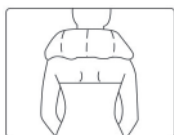

Neck

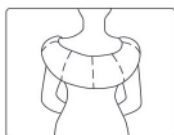

Back



Waist



Abdomen



Thighs

Warranty Information

Training equipment from TAURUS® is subject to strict quality control. However, if a fitness equipment purchased from us does not work perfectly, we take it very seriously and ask you to contact our customer service as indicated. We are happy to help you by phone via our service hotline.

Warranty period

The warranty periods, shown on our web shop, begin on delivery of the fitness equipment. The respective warranty periods for your equipment can be found on its product website.

Please refer to the website of your TAURUS® device for the exact warranty conditions.

Care, Maintenance and Storage

WARNING

The storage location should be chosen so that improper use by third parties or children is prevented.


1. Clean the device with a dry soft cloth.
2. Store the device in the supplied carton in a dry environment.

Disposal

At the end of its operational life, this equipment cannot be disposed of in normal household waste. Instead, it must be disposed of via a recycling centre. Further information can be obtained from your local authority's recycling service.

To ensure proper disposal of the massage device, the li-on battery must be removed and disposed of separately.

1. Open the zip along the upper side.
2. Disconnect the cable connection and remove the battery.

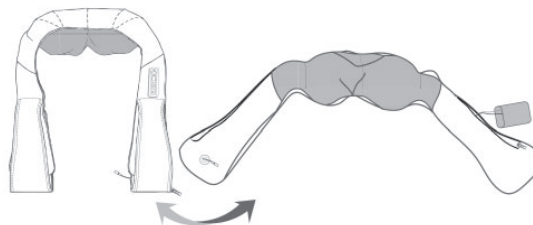

Contact

DE TECHNIK

+49 4621 4210-900

technik@fitshop.de

SERVICE

0800 20 20277 (kostenlos)

info@fitshop.de

UK CUSTOMER SUPPORT

00800 2020 2772

+494621 4210 944

info@fitshop.co.uk

FR TECHNIQUE & SERVICE

+33 (0) 189 530984

+49 4621 42 10 933

info@fitshop.fr

DK TEKNIK OG SERVICE

80 90 16 50

+49 4621 4210-945

info@fitshop.dk

ES TECNOLOGÍA Y SERVICIOS

911 238 029

info@fitshop.es

BE TECHNIQUE & SERVICE

02 732 46 77

+49 4621 4210 933

info@fitshop.be

NL TECHNISCHE DIENST & SERVICE

+31 172 619961

service@fitshop.nl

AT TECHNIK & SERVICE

0800 20 20277 (kostenlos)

+49 4621 42 10-0

info@fitshop.at

CH TECHNIK & SERVICE

0800 202 027

+49 4621 42 10-0

info@fitshop.ch

PL DZIAŁ TECHNICZNY I SERWIS

22 307 43 21

+49 4621 42 10-948

info@fitshop.pl

INT TECHNICAL SUPPORT & SERVICE

+49 4621 4210-944

service-int@fitshop.de